# PLUG & PRESS® DISPENSER

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Instructions for use

Mode d'emploi

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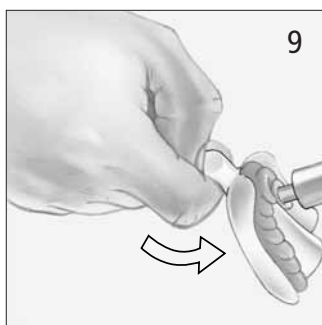
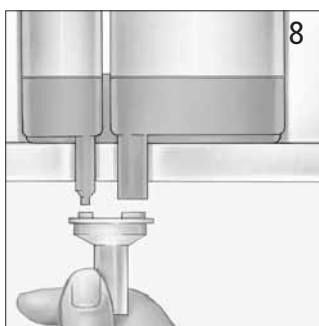
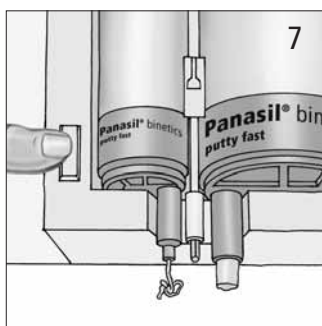
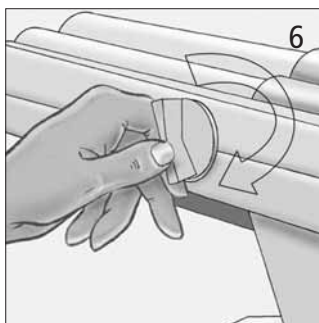
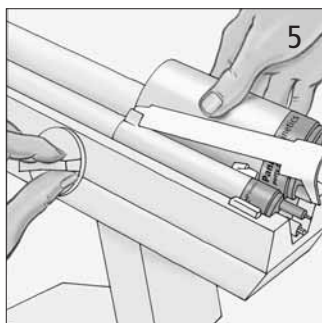
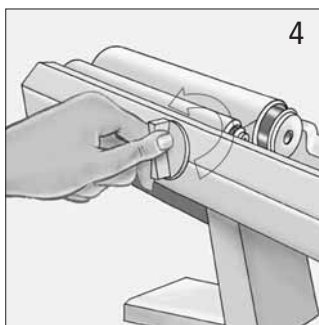
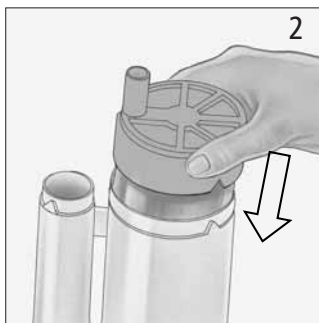
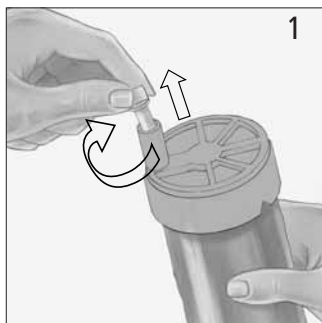
D

Deutsch

GB

English

F

Française

I

Italiano

E

Español

P

Português

NL

Nederlands

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# PLUG & PRESS® DISPENSER

## Automatic dispensing and mixing unit

Instructions for use

For professional use only. Caution: Federal (U.S.A.) law restricts this device to sale by or on the order of a dentist.

## Contents

	Page
1. Introduction	7
2. Package contents	7
3. Technical data	7
4. Safety instructions	7 - 8
5. Installation	8
6. Operation	8 - 9
6.1. Using foil bags	8 - 9
6.2. Using prefilled cartridges	9
7. Changing cartridges	9
8. Disposal of used materials	9
9. Care and maintenance	9
10. Storage / Transport	9
11. Faults / Troubleshooting	10
Terms of guarantee	31

## 1. Introduction

The Plug & Press® Dispenser is an automatic dispensing and mixing unit for use with impression materials in foil bags or cartridges in 1:5 systems. Suitable cartridge bodies are available for inserting the foil bag into the unit. The two components are mixed optimally using a dynamic mixer, which is also available as an accessory.

## 2. Package contents

Dispensing and mixing unit  
Mains adaptor  
Power cable  
Instructions for use  
Guarantee

## 3. Technical data

Drive: Low-voltage direct current motor  
Operating mode: Short-term operation  
Housing: All housing components are made from tough, unbreakable plastic  
Cartridges: Suitable for prefilled foil bags in cartridges or prefilled cartridges for 1:5 systems  
Diameter: 60 x 130 mm and 26.8 x 130 mm  
Mains voltage: External mains adaptor:  
    input 100-240V / 1.5A / AC / 50Hz-60Hz  
    output 48V / 1,0A / DC  
Unit connection:  
    48V / 1,0A / DC  
Frequency: 50 Hz – 60 Hz (external mains adaptor)  
Power consumption: 1.0 A (external mains adaptor)  
Hydraulic system: Operating pressure limited to max. 35 bar (by safety valve)  
Ambient temperature: 18 °C – 40 °C/64 °F – 104 °F  
Relative humidity: 20% – 80%  
Protective system: IP 4X  
Dimensions: L 370 mm x B 180 mm x H 290 mm  
Weight: 9.1 kg

## 4. Safety instructions

Read these instructions for use carefully before operating the unit.

1. Adhere strictly to the following instructions when operating the unit.  
We do not accept any liability for damage to the unit or a third person caused by incorrect use of this unit.
2. Before operating the unit, ensure that the mains voltage given on the adaptor corresponds with that

of the existing mains voltage. Connection to the incorrect mains voltage can damage the unit and prevent the unit operating properly.

3. Use the unit always by room temperatures (18 – 24 °C/64 – 76 °F).  
If the unit is taken from a cold area into a warm one, allow a period of time before switching on the unit. Condensation can be dangerous and cause a short circuit. This is why the unit should only be operated once it is at the ambient temperature.
4. Always ensure that the unit is connected to a socket with the same voltage(s) as displayed on the adaptor before operating the unit.
5. Position the unit on a nonslip, flat, level surface.
6. Always use suitable material in foil bags in the system cartridge body or prefilled cartridges.
7. Always ensure that foil bags or cartridges have been properly activated before operating the unit.
8. Do not insert any objects into the unit unless replacing components in accordance with these instructions for use.
9. Always use original components when replacing defect plunger discs according to these instructions for use. The manufacturer does not accept any liability for damage caused by the use of plunger discs from another manufacturer.
10. If the unit does not operate after following the above instructions or is visibly damaged, it should be assumed that the safety of the unit has been impaired. In this case the unit should be immediately taken out of service and clearly marked to ensure that it is not inadvertently operated again by a third party.
11. Ensure that solvents, flammable liquids and intense sources of heat are kept away from the unit. These could cause damage to the plastic housing or components inside the unit.
12. When cleaning the unit, do not use cleaning agents, as these may penetrate the housing and trigger an electrical short circuit or a hazardous malfunction.
13. The unit should only be opened and repaired by trained personnel.
14. **Do not start the unit if there is hardened material in the dynamic mixer, as this could damage the unit.**

## 5. Installation

The unit is designed for use on a nonslip, flat, level surface.

If the unit is moved from a very cold environment (< 5 °C) into a warm environment, wait approx. 2 hours before switching on the unit. The hydraulic oil in the unit should reach room temperature to ensure the unit operates properly.

Before operating the unit, ensure that the mains voltage given on the adaptor corresponds with that of the existing mains voltage.

The unit is ready for operation after connection to the mains supply.

**An electronic switch prevents the unit being switched on before the plungers are inserted into the cartridges.**

**6. Operation** (Please take note of the corresponding drawings on the inside cover.)

### 6.1. Using foil bags

Before using for the first time, remove the white safety pin from the activating head of the large foil bag by turning it in the direction of the arrow and pulling it out. (Figure 1)

Insert the large foil bag with the integrated activating head into the cartridge body. Ensure that the notches on the activating head and cartridge body are aligned. (Figure 2)

Press the activating head onto the cartridge body until it is fully closed. The foil bag is automatically pierced by a pin when the activating head is pressed in position. (Figure 3)

Follow the same procedure as in step 2 and step 3 with the smaller foil bag.

Move the plungers right to the top by turning the control knob. (Figure 4)

Insert the cartridge body into the unit. (Figure 5)

Turn the control knob to move the plungers into the cartridge body and continue turning until the plungers come into contact with the foil bags and the seal plungers of the cartridge. Once there is contact, press one of the two start keys (at the front on top of the unit) to dispense the material. (Figure 6)

Before fitting the dynamic mixer, extrude the material until it exits uniformly. (Figure 7)

Do not fit the dynamic mixer and close the locking clip until this has been completed. (Figure 8)

Fill the tray with the required amount of material. Hold the impression tray at an angle and press lightly against it. Leave the mixer in the material when dispensing. Leave the filled mixer on the cartridge body as a seal.

(Figure 9)

Before using again, release the locking clip to remove the used mixer and check that the outlets in the activating heads are not blocked. Fit a new mixer, close the locking clip and continue as usual.

Simply remove the empty foil bags and activating heads from the cartridge body and dispose of them. The cartridge body is reused. (Figure 10)

### 6.2. Using prefilled cartridges

Proceed as described in 6.1 for stages 4 – 9. The empty cartridge should be disposed of and should not be reused.

## 7. Changing cartridges

Turn the control knob on the side of the unit to move the plungers up to the start position. The cartridge body can only be removed from the unit when the plungers are at the top.

To insert the new cartridge body proceed as described in section 6. Operation.

## 8. Disposal of used materials

Dispose of used dynamic mixers and empty foil bags or cartridges with domestic waste.

## 9. Care and maintenance

With normal use the unit is virtually maintenance free.

### Care:

All parts should be cleaned with a soft cloth, which can be slightly dampened if required. No moisture should get into the unit, however, as this could cause a short circuit. Never use solvents and/or abrasive cleaning agents, as these damage the plastic components.

### Changing the plunger discs:

The plunger discs are parts subject to wear and tear and have to be changed from time to time.

Always remove the plug from the mains socket before changing the plunger discs.

For safety reasons always use original accessories.

Loosen the screw in the centre of the disc with an Allen key.

Remove the disc and screw a new plunger disc tightly in position.

Fit the plunger discs with the convex side facing the foil bag or the seal plungers of the cartridge.

### Note:

Do not open the unit housing. There are no parts in the housing that require maintenance or repair by the operator.

The manufacturer only accepts responsibility for the safety, reliability and efficient operation of the unit if it is used in accordance with these instructions for use. Cartridges are parts subject to wear and tear and are not covered by the guarantee.

## 10. Storage/Transport

Store/ transport at -20 °C to +60 °C/-4 °F to 140 °F. Maximum relative humidity 80%.

Do not throw away the unit in the household garbage.

The information provided for Kettenbach products is based on comprehensive research and experience in application technology.

Results are furnished to the best of our knowledge, subject to technical changes within the framework of product development.

However, users must comply with and consider all recommendations and information in connection with any use.

Updated: February 2010

## 11. Faults/Troubleshooting

Fault	Possible cause	Remedy
Unit does not start.	Unit not connected to electricity supply.	Check that the electricity cable is correctly connected and the plug is in the socket.
	Unit is not level.	Set the unit on a flat, level surface.
	Plungers have not been slid into the cartridge body.	Slide the plungers in using the hand wheel.
Unit starts, but no material is dispensed.	Foil bag is empty.	Insert a new foil bag into the cartridge.
	Cartridge is empty.	Insert a new cartridge.
	Foil bag is not open or not sufficiently open.	Open the bag properly (see "Operation" steps 1, 2 and 3).
	Hardened material in an old dynamic mixer.	Fit new dynamic mixer.
	Plunger does not reach foil bag, as the foil bag is half empty.	Move the plunger down manually (see "Operation" step 6).
	Plunger jammed.	Check ease of movement by turning the control knob. Change the foil bag if necessary.
	Foil bag inserted incorrectly.	Correct the position or re-insert the foil bag correctly (see "Operation" steps 1, 2 and 3).
Material is dispensed, but the dynamic mixer does not rotate.	Drive shaft does not engage in the dynamic mixer.	Attach the dynamic mixer properly.
Cartridge cannot be inserted into the unit.	Plungers are not right at the top.	Move the plungers right to the top by turning the control knob anti-clockwise. Check the plunger discs for damage.
	Activating head in the wrong position or foil bag not properly activated.	Check the fit of the activating head. Correctly activate the foil bag if necessary (see "Operation" steps 1, 2 and 3).
Cartridge cannot be removed from the unit.	Plungers are not right at the top.	Move the plungers right to the top by turning the control knob anti-clockwise. Check the plunger discs for damage.
Unit does not switch off.	Start key jammed.	Press the start key again, remove the mains plug if necessary and remove the obstruction from the start key.
Plunger difficult to retract.	Foil bag jammed between the cartridge and the plungers.	Move the plungers away from the obstruction, avoid damaging the plunger discs. Use a new foil bag of impression material. If necessary change the plunger discs and use a new cartridge body.
	Unit has no power supply.	Check that the electricity cable is correctly connected and the plug is in the socket.
	Unit is not level.	Set the unit on a flat, level surface.
Hydraulic liquid escapes from the unit.	Hydraulic system not properly sealed.	Do not operate the unit. Inform your dental dealer for help.

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Distributed in the U.S. by:  
**Kettenbach LP**  
7777 Center Avenue, Suite 280  
Huntington Beach, CA 92647  
USA

**KETTENBACH SNC**  
7 Place de la Gare  
57200 Sarreguemines  
Tél.: 08100 200 60  
Fax: 0800 800 977

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please take note of the inside cover

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## Terms of guarantee

GB

Kettenbach guarantee perfect workmanship and quality of this appliance for a period of 24 months - excluding parts subject to normal wear and tear. This guarantee is null and void if inspection proves that the damage is due to unauthorized interference, incorrect operation or non-observance of the operating instructions. This includes the use of materials not supplied by Kettenbach, provided they are not the cause of the defect resulting in the complaint. The purchaser is responsible for proving that such materials are not the cause of complaint. The guarantee is fulfilled by repairing or replacing the appliance free of charge. Further claims are excluded, especially compensation for additional damage, direct or indirect. Statutory guarantee claims are not affected hereby.

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